
# Mepilex<sup>®</sup> Border Lite

WITH SAFETAC<sup>®</sup> TECHNOLOGY

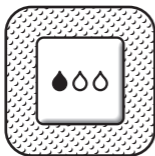STERILE EO

CE 2797

MÖLNLYCKE  
HEALTH CARE

Manufacturer

Mölnlycke Health Care AB

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# Mepilex® Border Lite

Self-adherent soft silicone foam dressing

## Product description

Mepilex Border Lite is a self-adherent, absorbent dressing that maintains a moist wound environment. The waterproof outer layer protects the wound from dirt and bacteria. The dressing has a Safetac® wound contact layer that is a unique adhesive technology. It minimises pain to patients and trauma to wound and skin at dressing removal.

Mepilex Border Lite consists of:

- a wound contact layer consisting of soft silicone adhesive (Safetac) and a film carrier
- a flexible absorbent pad in two layers: a foam and a non-woven spreading layer
- an outer film which is breathable but waterproof, providing a barrier to external contaminants

Dressing material content:

Silicone, polyurethane, polyacrylate, viscose, polyester and polyolefin.

## Indications for use

Mepilex Border Lite is designed for the management of a wide range of non/low exuding wounds, such as leg and foot ulcers, pressure ulcers, surgical wounds and traumatic wounds e.g. abrasions, blisters and skin tears.

Mepilex Border Lite can also be used as protection of compromised and/or fragile skin.

## Precautions

- Do not use on patient with known hypersensitivity to the dressing or its ingoing materials.
- Do not use together with oxidising agents such as hypochlorite solutions or hydrogen peroxide.
- If you see signs of clinical infection e.g. fever or the wound or surrounding skin becoming red, warm or swollen, consult a health care professional for appropriate treatment.

- Do not reuse. If reused performance of the product may deteriorate, cross contamination may occur.
- Sterile. Do not use if sterile barrier is damaged or opened prior to use. Do not re-sterilise.

## Instructions for use

Mepilex Border Lite may be used by lay persons under supervision of health care professionals.

1. Cleanse the wound according to clinical practice. Dry the surrounding skin thoroughly.
2. Select an appropriate dressing size/shape. The wound pad should cover the dry surrounding skin by at least 1-2cm.
3. Remove the first release film and apply the adherent side to the wound.
4. Remove the remaining release film(s) and smooth down the border on the skin. Do not stretch the dressing.

The dressing change interval may be several days. Change the dressing before it is fully saturated, at signs of leakage, or as indicated by clinical practice. Mepilex Border Lite can be used under compression bandaging, and in combination with gels.

## Special storage conditions and handling conditions

Mepilex Border Lite should be stored in dry conditions.

The foam may change colour to more yellow when exposed to light, air and/or heat. This has no influence on product properties.

Disposal should be handled according to local environmental procedures.

## Other information

If any serious incident has occurred in relation to the use of Mepilex Border Lite, it should be reported to Mölnlycke Health Care.

Mepilex® and Safetac® are registered trademarks of Mölnlycke Health Care AB.

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